



- Extensive Range
- 90+ Years of History
- Extended Aftercare Support
- UKAS Laboratory

Intellcold® Pharmacy Fridges

Free Standing RLDF1519

A professional 440 litre pharmacy refrigerator powered by our unique IntelliCold® touch screen controller, designed from the ground up for the safe storage of cold chain medicines as advised by Public Health England in Protocol for Ordering, Storing and Handling Vaccines.

It features a 1.6 metre secondary (wandering) temperature probe that can be placed anywhere inside the pharmacy fridge and it features a NEW touch screen controller which is completely flat so can be thoroughly cleaned and disinfected as there are no raised buttons to collect dirt and germs. The controller displays, alarms and records minimum and maximum temperatures for both temperature probes. Like a mobile phone touch screen it locks after 30 seconds to ensure that buttons can't be activated accidentally.

We can also provide a calibration certificate for your pharmacy refrigerator from our in house UKAS calibration laboratory No: 8898 as an optional extra.

We are confident that this pharmacy refrigerator will provide you with many years of fault free service. However, for extra peace of mind your RLDF1519 comes complete with a FREE UK 5 year parts and labour warranty.

Specifications

Product Code	Free Standing RLDF1519
Gross Capacity	440 litres
Style	Free Standing
UK Warranty	5 years
External Dimensions (mm) HxWxD	1875 x 600 x 700
Internal Dimensions (mm) HxWxD	1724 x 480 x 500
Temperature Range	2°C to 8°C
Fridge Set Point	5°C

Specifications (cont.)

Internal Material	Powder Coated Steel
External Material	Powder Coated Steel
Shelves	6
Net Weight	93kg
Door Style	Solid
Lock Type	Key Lock
Digital Temperature Display	Yes
Digital Temperature Control	Yes
High/Low Temperature Alarm	Yes
Min/Max Temperature Recording	Yes
Integral Temp Data Logging	Yes
Automatic Temperature & Fridge Operation Data	Yes
Number of Temperature Probes	2
Internal Light	Yes
Self Closing Door	Yes
Door Open Alarm	Yes
Mains Failure Alarm	Yes
Alarm Battery Backup	Yes
Auto Defrost	Yes
Fan Assisted Cooling	Yes
Remote Alarm Contacts	Yes
Access Porthole	Yes
Access Porthole Size	12mm
Power Consumption Kwh/24hrs	2.55
Ambient Operating Temperature	16°C to 32°C

Features

- Fitted with an IntelliCold® touch screen controller
- Controller and alarm battery back up
- Power failure alarm
- Low temperature alarm
- Off cycle auto defrost
- Access port
- Ammonia free
- Self closing door
- Secondary (wandering) temperature probe
- Remote alarm contacts
- Compliant with all leading professional guidelines
- 24/7 automatic temperature data logging which is accessible via transferring data to your computer
- Door open alarm
- High temperature alarm
- Minimum/maximum temperature recording and display
- Fan assisted cooling
- Environmentally friendly refrigerant
- Door lock with two keys
- Real time temperature display for two temperature probes
- Internal light
- FREE five year labour warranty

Models

	Model	Gross Capacity	Door Style	Lock Type
	Bench Top RLDF0119 Diglock	36 litres	Solid	Digital Lock
	Bench Top RLDG0219 Diglock	66 litres	Glass	Digital Lock
	Bench Top RLDF0224	66 Litres	Solid	Key Lock
	Bench Top RLDF0224DIGLOCK	66 Litres	Solid	Digital Lock
	Bench Top RLDF0124	36 litres	Solid	Key lock
	Bench Top RLDG0124	36 litres	Glass	Key Lock
	Bench Top RLDG0224	66 Litres	Glass	Key Lock
	Bench Top RLDG0124DIGLOCK	36 litres	Glass	Digital Lock
	Free Standing RLDF1519	440 litres	Solid	Key Lock

Accessories



Micro SD Card




Divided Drawer for
RLDF/G1019,
RLDF/G1519,
RLBM1019 and
RLBM1519



Data Logger
RMBL2019-M1



IntelliCold® Remote
Alarm

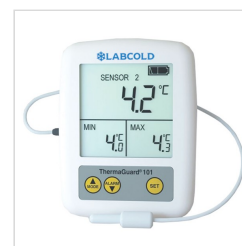

Digital Thermometer
and Alarm
RLAA5003-1