



LABCOLD™

medical & scientific refrigeration




RDBB2320MD Dual Compressor Blood Bank

Part Number	RDBB2320MD
Capacity (Litres)	485
Capacity (250ml bags)	425 (approx)
Style	Freestanding
Exterior Dimensions (HxWxD)	1990 x 716 x 855mm
Interior Dimensions (HxWxD)	1410 x 470 x 565mm
Weight	183kg
Energy Consumption	5.5 kWh/24h
Operating Temperature	+4°C +/- 2°C
Interior Construction	Stainless Steel
Exterior Construction	Powder Coated Steel
Door Style	Solid
Door Hinging	Right
Dual Refrigeration	✓
Open Door Alarm	✓
Load Temperature Alarm	✓
Air Temperature Alarm	✓
Mains Failure Alarm	✓
Air Temperature Display	✓
Load Temperature Display	✓
7 Day Chart Recorder	✓
Lock	Digital Keypad
Battery Back Up	✓
Drawers	6
Remote alarm contacts	✓
Access Port	✓
Auto Defrost	✓
Internal Light	✓
Eco Friendly Refrigerant	✓
Certified Medical Device in Accordance with UK MDR 2002	✓

Fitted with a Dual Compressor to ensure greater temperature control and peace of mind in the event of failure of one of the systems. Each complete refrigeration system provides 100% performance, sufficient to cool down and maintain temperatures correctly without help from the second compressor.



Hygienic touch screen controller

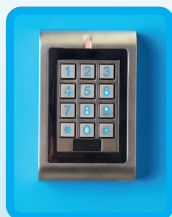

Digital keypad lock

An asymmetric change over ensures even distribution of workload, resulting in even wear of the compressors and also allows for one to defrost while the other is still running. The two independent systems provide back up should an issue develop plus there is no need to move contents immediately if a problem is identified.

Designed and manufactured in the UK in our Basingstoke factory, this Blood Bank refrigerator is a certified Medical Device in accordance with UK MDR 2002 and is compliant with BS4376:Part 1:1991.

Labcold blood banks are tested for 72 hours and calibrated by our UKAS accredited calibration laboratory #8898 before they are dispatched and all come with a FREE 2 year parts and labour warranty.

01256 705 570

www.labcold.com

