


## RDBG2068MD Compact Dual Compressor Blood Bank



Ideal for A&E, Maternity, Community Hospitals, and Hospices etc, fitted with a Dual Compressor to ensure greater temperature control and peace of mind in the event of failure of one of the systems. Each complete refrigeration system provides 100% performance, sufficient to cool down and maintain temperatures correctly without help from the second compressor.



Hygienic touch screen controller

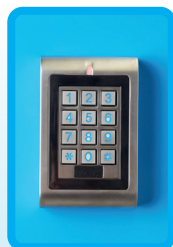

Digital keypad lock

An asymmetric change over ensures even distribution of workload, resulting in even wear of the compressors and also allows for one to defrost while the other is still running. The two independent systems provide back up should an issue develop plus there is no need to move contents immediately if a problem is identified.

Designed and manufactured in the UK in our Basingstoke factory, this Blood Bank refrigerator is a certified Medical Device in accordance with UK MDR 2002 and is compliant with BS4376:Part 1:1991.

Labcold blood banks are tested for 72 hours and calibrated in our UKAS calibration laboratory No:8898 before they are dispatched and all come with a FREE 2 year parts and labour warranty.

| Part Number   | RDBG2068MD          |
|---|---------------------|
| Capacity (Litres)                                       | 118                 |
| Capacity (250ml bags)                                   | 68 (approx)         |
| Style   | Freestanding        |
| Exterior Dimensions (HxWxD)                             | 1045 x 655 x 670mm  |
| Interior Dimensions (HxWxD)                             | 445 x 560 x 475mm   |
| Weight  | 85kg                |
| Energy Consumption                                      | 2.4 kWh/24h         |
| Operating Temperature                                   | +4°C +/- 2°C        |
| Interior Construction                                   | Stainless Steel     |
| Exterior Construction                                   | Powder Coated Steel |
| Door Style  | Glass               |
| Door Hinging  | Right               |
| Dual Refrigeration                                      | ✓                   |
| Open Door alarm   | ✓                   |
| Load Temperature Alarm                                  | ✓                   |
| Air Temperature Alarm                                   | ✓                   |
| Mains Failure Alarm                                     | ✓                   |
| Air Temperature Display                                 | ✓                   |
| Load Temperature Display                                | ✓                   |
| 7 Day Chart Recorder                                    | ✓                   |
| Lock  | Digital Keypad      |
| Battery Back up (Display)                               | ✓                   |
| Drawers   | 2                   |
| Remote alarm contacts                                   | ✓                   |
| Access Port   | ✓                   |
| Auto Defrost  | ✓                   |
| Internal Light  | ✓                   |
| Eco Friendly Refrigerant                                | ✓                   |
| Certified Medical Device in Accordance with UK MDR 2002 | ✓                   |